

NG

Технические характеристики

По вопросам продаж и поддержки обращайтесь:

Алматы (7273)495-231
Архангельск (8182)63-90-72
Астрахань (8512)99-46-04
Барнаул (3852)73-04-60
Белгород (4722)40-23-64
Брянск (4832)59-03-52
Владивосток (423)249-28-31
Волгоград (844)278-03-48
Вологда (8172)26-41-59
Воронеж (473)204-51-73
Екатеринбург (343)384-55-89
Иваново (4932)77-34-06
Ижевск (3412)26-03-58
Иркутск (395)279-98-46
Россия (495)268-04-70

Казань (843)206-01-48
Калининград (4012)72-03-81
Калуга (4842)92-23-67
Кемерово (3842)65-04-62
Киров (8332)68-02-04
Краснодар (861)203-40-90
Красноярск (391)204-63-61
Курск (4712)77-13-04
Липецк (4742)52-20-81
Магнитогорск (3519)55-03-13
Москва (495)268-04-70
Мурманск (8152)59-64-93
Набережные Челны (8552)20-53-41
Нижний Новгород (831)429-08-12
Киргизия (996)312-96-26-47

Новокузнецк (3843)20-46-81
Новосибирск (383)227-86-73
Омск (3812)21-46-40
Орел (4862)44-53-42
Оренбург (3532)37-68-04
Пенза (8412)22-31-16
Пермь (342)205-81-47
Ростов-на-Дону (863)308-18-15
Рязань (4912)46-61-64
Самара (846)206-03-16
Санкт-Петербург (812)309-46-40
Саратов (845)249-38-78
Севастополь (8692)22-31-93
Симферополь (3652)67-13-56
Казахстан (7172)727-132

Смоленск (4812)29-41-54
Сочи (862)225-72-31
Ставрополь (8652)20-65-13
Сургут (3462)77-98-35
Тверь (4822)63-31-35
Томск (3822)98-41-53
Тула (4872)74-02-29
Тюмень (3452)66-21-18
Ульяновск (8422)24-23-59
Уфа (347)229-48-12
Хабаровск (4212)92-98-04
Челябинск (351)202-03-61
Череповец (8202)49-02-64
Ярославль (4852)69-52-93

Histology

Bath



Paraffin Section Mounting Bath

Comprising aluminium housing with PTFE coated matt black interior, illuminated on/off rocker switch, rotary control dial, energy regulator and resettable thermal fuse providing overtemperature protection. Without lid.

Catalogue No.	HG300-45	
Model	MH8517	
Control	Thermostatic	
Range	°C	Ambient to 70, $\pm 0.5^*$
Bath volume	litres	2.5
depth	mm	80
diameter	mm	220
Rating	W	240

*When filled with 2 litres of water, and after 120 minutes.

As described. Overall 330 x 140mm diameter x height. Weight 2.7kg. For 230V 50/60Hz single phase supplies.

HG300-45 Ref. MH8517

HG300-55 Lid, black PP, for MH8516



HG300



HG310 in use

Drying Bench



Slide Drying Bench

With DHSS ESCHLE approval. With built in energy regulator and indicator lamps. Maximum temperature 100°C. Accepts up to fifty 76 x 25mm slides on removable drying bars. Overall 435 x 178 x 75mm W x D x H. For 220-240V a.c. single phase supplies, 150W.

HG310-10 Ref. MH6616

Slide rack ends

For suspending microscope slides over laboratory sinks using glass rod rails. Each end has two 7.5mm diameter holes with screw clamps to accommodate the glass rods, and two levelling screw feet. Supplied in pairs.

ML754-10 Staining rack ends

Other Histology Items

Staining Troughs/Jars, Microscope Slides – see *Microscopy section*.

Forceps/Scalpels – see *Dissecting section*.

Wax Dispenser



Paraffin Wax Dispenser

With 7.5 litre capacity, stainless steel tank and rapid heating system for melting up to 6kg of pelletised wax. Operating range: ambient to +70°C. The non-dripping, lever operated tap has three modes, off, on-demand dispense and continuous flow and is separately heated to avoid blockage. The tank outlet also features a 0.5mm mesh filter to prevent coarse particles fouling the tap. Digital setting and control with $\pm 1^\circ\text{C}$ accuracy, safety circuit breaker and temperature cut-out at 120°C. Overall 455 x 181 x 500mm H x W x D, weight 7kg. For 220-240V 50/60Hz a.c. single phase supplies, 480W.

HG320-20 Ref. MH8524



HG320



ML754 in use

Histology

Tissue-Wax

Tissue Embedding Wax, Tissue-Wax

A highly purified, filtered, white paraffin wax blended with synthetic polymers and suitable for routine embedding and sectioning applications. Offers very little shrinkage on cooling and high performance when sectioning, ease of cutting and ribbon continuity. Melting point 56-58°C, congealing point 54-57°C. Supplied in pelletised form within packs of 1kg.

HG385-15 Tissue-Wax

Tissue Embedding Wax, Tissue-Wax*

Similar to HG385-15 but with a small amount of Dimethyl Sulphoxide added to enhance tissue penetration, giving a shorter infiltration time. Also recommended for larger or harder tissue samples as it reduces brittleness of sections, enabling more flexible ribbons. Supplied in pelletised form within packs of 1kg.

HG385-25 Tissue-Wax*

Larger pack sizes are available - *details on request.*

Embedding Cassettes

Embedding Cassettes, Standard

Solvent-resistant polyoxymethylene (POM). Made to fit most commonly-used microtome chucks. Circular base perforations enhance adhesion of paraffin wax blocks to the cassette. Without cover, but with moulded lugs to accept accessory stainless steel covers. A 45° angled labelling area provides a comfortable writing surface and is suitable for most models of cassette labeller. Supplied in packs of 1000 (2 x bags of 500), in colours as indicated.

HG420-02 Light blue
HG420-04 Yellow
HG420-06 Light green
HG420-08 Grey

HG420-10 White
HG420-12 Pink
HG420-14 Beige
HG420-16 Orange

Embedding Cassettes, Universal

Generally as HG420 but with slotted perforations on base and fitted, hinged lid. Supplied in packs of 500, in colours as indicated.

HG424-03 Light blue
HG424-05 Yellow
HG424-07 Light green
HG424-09 Grey

HG424-11 White
HG424-13 Pink
HG424-15 Beige
HG424-17 Orange

Biopsy Cassettes


Biopsy Cassettes

Polyoxymethylene (POM). With finely perforated (1mm x 1mm aperture) base and fitted hinged lid. Supplied in packs of 500, in colours as indicated.

HG429-04 Light blue
HG429-07 Yellow
HG429-10 Light green
HG429-13 Grey

HG429-16 White
HG429-19 Pink
HG429-22 Beige
HG429-25 Orange

Marker pens for biopsy cassettes - see LA620-series.

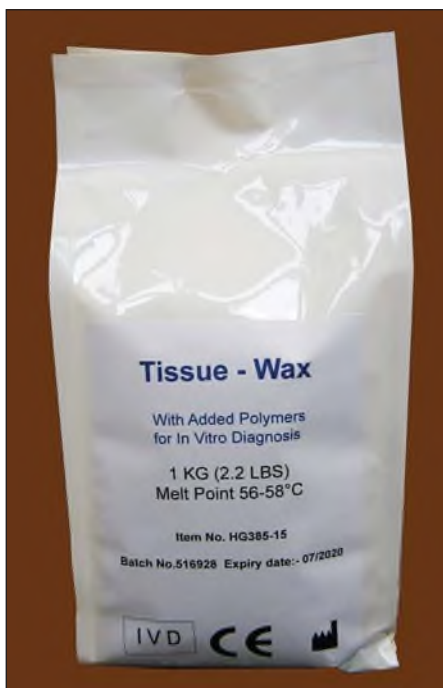

HG385-15, HG385-25 is similar



HG420



HG424



HG429

По вопросам продаж и поддержки обращайтесь:

Алматы (7273)495-231	Казань (843)206-01-48	Новокузнецк (3843)20-46-81	Смоленск (4812)29-41-54
Архангельск (8182)63-90-72	Калининград (4012)72-03-81	Новосибирск (383)227-86-73	Сочи (862)225-72-31
Астрахань (8512)99-46-04	Калуга (4842)92-23-67	Омск (3812)21-46-40	Ставрополь (8652)20-65-13
Барнаул (3852)73-04-60	Кемерово (3842)65-04-62	Орел (4862)44-53-42	Сургут (3462)77-98-35
Белгород (4722)40-23-64	Киров (8332)68-02-04	Оренбург (3532)37-68-04	Тверь (4822)63-31-35
Брянск (4832)59-03-52	Краснодар (861)203-40-90	Пенза (8412)22-31-16	Томск (3822)98-41-53
Владивосток (423)249-28-31	Красноярск (391)204-63-61	Пермь (342)205-81-47	Тула (4872)74-02-29
Волгоград (844)278-03-48	Курск (4712)77-13-04	Ростов-на-Дону (863)308-18-15	Тюмень (3452)66-21-18
Вологда (8172)26-41-59	Липецк (4742)52-20-81	Рязань (4912)46-61-64	Ульяновск (8422)24-23-59
Воронеж (473)204-51-73	Магнитогорск (3519)55-03-13	Самара (846)206-03-16	Уфа (347)229-48-12
Екатеринбург (343)384-55-89	Москва (495)268-04-70	Санкт-Петербург (812)309-46-40	Хабаровск (4212)92-98-04
Иваново (4932)77-34-06	Мурманск (8152)59-64-93	Саратов (845)249-38-78	Челябинск (351)202-03-61
Ижевск (3412)26-03-58	Набережные Челны (8552)20-53-41	Севастополь (8692)22-31-93	Череповец (8202)49-02-64
Иркутск (395)279-98-46	Нижний Новгород (831)429-08-12	Симферополь (3652)67-13-56	Ярославль (4852)69-52-93
Россия (495)268-04-70	Киргизия (996)312-96-26-47	Казахстан (7172)727-132	